

bioloren[®]

2025/2026

PRODUCT CATALOG

metal free dental solutions



AVANT
Fiber Posts

REAMERS
Endodontic

AVANT
Special Posts

InFibra
Ribbon

URC
Cements

TRILOR
Discs, Arches
and Blocks

PPE
Personal
Protection
Equipment

MADE IN ITALY



OUR COMPANY

BILOREN is an innovative Italian company specializing in "Metal free" dental systems.

The company was founded in 1998 and has since focused on the manufacture of cutting-edge dental products.

Over the years the company has earned a considerable reputation as a producer of fiber posts sold in more than 40 countries.

The quality system of BIOLOREN is certified UNI CEI EN ISO 13485. All the products by BIOLOREN are certified CE and many of them also received the prestigious American FDA.

BIOLOREN uses high quality materials and cooperates with important universities both in Italy and abroad.

The products are designed and manufactured by BIOLOREN which owns patents for many of them. In particular, the research team specializes in the development of fiberglass, carbon and polyethylene products.

BILOREN provides metal free solutions to meet the needs of all dentists and dental technicians. The company distributes worldwide through a network of distributors.



Certificate of Registration

QUALITY MANAGEMENT SYSTEM - ISO 13485:2016

This is to certify that:

BIOLORENS Srl

Via Vela, 50

Saronno (VA)

21042

Italy

Has the Capability to:

FM 701537

and operates a Quality Management System which complies with the requirements of ISO 13485:2016 for the following scope:

Design, manufacture and distribution of non-sterile prosts, bars, blocks, ribbon for dentistry and related metal-free materials for dental prosthetics.

Previous certificate expires on 2022-03-26

Recertification audit ended 2022-03-21



For and on behalf of BSI:


Graeme Tunbridge, Senior Vice President Medical Devices

original Registration Date: 2015-04-16

Latest Revision Date: 2022-03-25

Effective Date: 2022-03-25

Expiry Date: 2025-03-26

Page: 1 of 1

..making excellence a habit™

This certificate remains the property of BSI and shall not be loaned, resold, rented, given, or disposed of in any manner without the written authority of BSI. For more information please visit www.bsigroup.com/bsi/certificates

To be used in conjunction with the scope shown on the attached schedule.

Approved by the BSI Group of Companies

Approved by the BSI Group of Companies Ltd., (2020) BSI House, 389, Market Road, Singapore, Tel: 65 6349 0001 Fax: 65 6349 0002

			
CERTIFICATO CE - GARANZIA DI QUALITÀ DELLA PRODUZIONE EC CERTIFICATE - PRODUCTION QUALITY ASSURANCE			
APPROVAZIONE DEL SISTEMA DI QUALITÀ ATTUATO DA APPROVAL OF THE QUALITY SYSTEM OPERATED BY			
BIOLOREN SRL			
IT - 21047 SARONNO (VA) - VIA ALESSANDRO VOLTA, 88			
BITI / BITES			
IT - 21047 SARONNO (VA) - VIA ALESSANDRO VOLTA, 88			
PER I SEGUENTI DISPOSITIVI O GRUPPI DI DISPOSITIVI / FOR THE FOLLOWING DEVICES OR GROUPS OF DEVICES <i>Decks, benches, tables and arches in fibre and resin for use outdoors.</i>			
<i>Decks, benches, tables and arches in fibre and resin for dental use</i>			
Certiquality S.r.l., Organismo Notificato n° 5546, certifica che il sistema di qualità Certiquality S.r.l., Notified Body n°5546, certifies that the quality system			
è conforme ai requisiti della Direttiva 93/42/CEE, Allegato is in compliance with the requirements of Directive 93/42/EEC, Annex		V	
RAPPORTO DI AUDIT N° AUDIT REPORT NO.	151682	CERTIFICATO N° CERTIFICATE N°	151682
Il presente certificato è emesso in duplice copia. Il originale, in possesso del richiedente, è il documento che ha valore dimostrativo di idoneità qualità e che può essere esibito a richiesta. Il duplicato, in possesso del notificatore, è il documento che attesta la validazione del sistema qualità e che può essere esibito a richiesta. Il presente certificato è emesso in duplice copia. Il originale, in possesso del richiedente, è il documento che ha valore dimostrativo di idoneità qualità e che può essere esibito a richiesta. Il duplicato, in possesso del notificatore, è il documento che attesta la validazione del sistema qualità e che può essere esibito a richiesta.			
La presente certificazione è valida per un periodo di validità di 3 anni, a partire dalla data di emissione del certificato. Il richiedente deve sottoporre il proprio sistema qualità a una revisione di idoneità qualità da parte del notificatore prima della scadenza della validità del certificato. Il notificatore ha il diritto di effettuare verifiche di idoneità qualità a sorpresa. Il richiedente deve sottoporre il proprio sistema qualità a una revisione di idoneità qualità da parte del notificatore prima della scadenza della validità del certificato. Il notificatore ha il diritto di effettuare verifiche di idoneità qualità a sorpresa. Il richiedente deve sottoporre il proprio sistema qualità a una revisione di idoneità qualità da parte del notificatore prima della scadenza della validità del certificato. Il notificatore ha il diritto di effettuare verifiche di idoneità qualità a sorpresa.			
PRIMA EMISSIONE FIRST ISSUE	23/06/2000		
REVISIONI SUCCESSIVE SUBSEQUENT ISSUES	23/06/2000		
DATA DI RINNOVO RENEWAL DATE	26/06/2004	R. PRESIDENTE - CESARE PUGGIONI	

Form 1-05-01-0001

CERTIFICAZIONE N° 1 EMISSIONE DI CERTIFICAZIONE DELLA QUALITÀ
 Via S. Gerardo, 2 - 20123 Milano - Tel. 02/58311111 - Fax 02/58311112 - certiquality@certiquality.it - www.certiquality.it

[illegible][illegible]

AVANT® FIBER POSTS

CODE Price	PRODUCT	SIZE	DESCRIPTION	TECNICAL INFORMATION
WPCY 08 WPCY 10 WPCY 12 WPCY 14 WPCY 17 WPCY 40	Fiber Post Cylindrical Glass 	ø 0,8 ø 1,0 ø 1,2 ø 1,4 ø 1,7 ø 0,8 ø 1,0 ø 1,2 ø 1,4	Fiberglass Cylindrical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. ø from 0,8 to 1,7 mm. Parallel shape.
TPCY 08 TPCY 10 TPCY 12 TPCY 14 TPCY 17 TPCY 40	Fiber Post Cylindrical Glass Translucent 	ø 0,8 ø 1,0 ø 1,2 ø 1,4 ø 1,7 ø 0,8 ø 1,0 ø 1,2 ø 1,4	Fiberglass Translucent Cylindrical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. ø from 0,8 to 1,7 mm. Parallel shape.
CPCY 08 CPCY 10 CPCY 12 CPCY 14 CPCY 17 CPCY 40	Fiber Post Cylindrical Carbon 	ø 0,8 ø 1,0 ø 1,2 ø 1,4 ø 1,7 ø 0,8 ø 1,0 ø 1,2 ø 1,4	Carbon Fiber Cylindrical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. ø from 0,8 to 1,7 mm. Parallel shape.
WPCO2 08 WPCO2 10 WPCO2 12 WPCO2 14 WPCO2 40	Fiber Post Conical Glass 2% 	Conicity 2% ø 0,8 - 1,25 ø 1,0 - 1,45 ø 1,2 - 1,65 ø 1,4 - 1,85	Fiberglass Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: ø from 0,8 to 1,4 mm. Upper: ø from 1,25 to 1,85 mm. Tapered shape 2%.
TPCO2 08 TPCO2 10 TPCO2 12 TPCO2 14 TPCO2 40	Fiber Post Conical Glass Translucent 2% 	Conicity 2% ø 0,8 - 1,25 ø 1,0 - 1,45 ø 1,2 - 1,65 ø 1,4 - 1,85	Fiberglass Translucent Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: ø from 0,8 to 1,4 mm. Upper: ø from 1,25 to 1,85 mm. Tapered shape 2%.
CPCO2 08 CPCO2 10 CPCO2 12 CPCO2 14 CPCO2 40	Fiber Post Conical Carbon 2% 	Conicity 2% ø 0,8 - 1,25 ø 1,0 - 1,45 ø 1,2 - 1,65 ø 1,4 - 1,85	Carbon Fiber Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: ø from 0,8 to 1,4 mm. Upper: ø from 1,25 to 1,85 mm. Tapered shape 2%.

AVANT® FIBER POSTS

CODE Price	PRODUCT	SIZE	DESCRIPTION	TECNICAL INFORMATION
WPCLRCO 08 WPCLRCO 10 WPCLRCO 12 WPCLRCO 14 WPCLRCO 40 WPCLRCO 40+4D	Fiber Post Conical Glass 2% Colour Code 	Conicity 2% $\varnothing$ 0,8 - 1,25 ○ $\varnothing$ 1,0 - 1,45 ● $\varnothing$ 1,2 - 1,65 ● $\varnothing$ 1,4 - 1,85 ●	Fiberglass Conical Post for endodontic reconstruction. Colour Code. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,8 to 1,4 mm. Upper: $\varnothing$ from 1,25 to 1,85 mm. Tapered shape 2%.
TPCLRCO 08 TPCLRCO 10 TPCLRCO 12 TPCLRCO 14 TPCLRCO 40	Fiber Post Conical Glass Translucent Colour Code 	Conicity 2% $\varnothing$ 0,8 - 1,25 ○ $\varnothing$ 1,0 - 1,45 ● $\varnothing$ 1,2 - 1,65 ● $\varnothing$ 1,4 - 1,85 ●	Fiberglass Translucent Conical Post for endodontic reconstruction. Colour Code. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,8 to 1,4 mm. Upper: $\varnothing$ from 1,25 to 1,85 mm. Tapered shape 2%.
WPCO2-4 12 WPCO2-4 14 WPCO2-4 16 WPCO2-4 18 WPCO2-4 40	Fiber Post Conical Glass 2-4% 	Conicity 2-4% $\varnothing$ 0,5 - 1,2 $\varnothing$ 0,7 - 1,4 $\varnothing$ 0,9 - 1,6 $\varnothing$ 1,1 - 1,8	Progressive conicity Fiberglass Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,2 to 1,8 mm. Tapered shape 2-4%.
TPCO2-4 12 TPCO2-4 14 TPCO2-4 16 TPCO2-4 18 TPCO2-4 40	Fiber Post Conical Glass 2-4% Translucent 	Conicity 2-4% $\varnothing$ 0,5 - 1,2 $\varnothing$ 0,7 - 1,4 $\varnothing$ 0,9 - 1,6 $\varnothing$ 1,1 - 1,8	Progressive conicity Fiberglass Translucent Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,2 to 1,8 mm. Tapered shape 2-4%.
WPCO4-6 15 WPCO4-6 17 WPCO4-6 19 WPCO4-6 21 WPCO4-6 40	Fiber Post Conical Glass 4-6% 	Conicity 4-6% $\varnothing$ 0,5 - 1,5 $\varnothing$ 0,7 - 1,7 $\varnothing$ 0,9 - 1,9 $\varnothing$ 1,1 - 2,1	Progressive conicity Fiberglass Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,5 to 2,1 mm. Tapered shape 4-6%.
TPCO4-6 15 TPCO4-6 17 TPCO4-6 19 TPCO4-6 21 TPCO4-6 40	Fiber Post Conical Glass 4-6% Translucent 	Conicity 4-6% $\varnothing$ 0,5 - 1,5 $\varnothing$ 0,7 - 1,7 $\varnothing$ 0,9 - 1,9 $\varnothing$ 1,1 - 2,1	Progressive conicity Fiberglass Translucent Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,5 to 2,1 mm. Tapered shape 4-6%.

AVANT® FIBER POSTS

CODE Price	PRODUCT	SIZE	DESCRIPTION	TECNICAL INFORMATION
TPCOGR 08 TPCOGR 10 TPCOGR 12 TPCOGR 14 0 TPCOGR 40 TPCOGR 40+4D	Fiber Post Conical Glass Translucent GRIPPER 	$\varnothing$ 0,8 - 1,25 $\varnothing$ 1,0 - 1,45 $\varnothing$ 1,2 - 1,65 $\varnothing$ 1,4 - 1,85	Fiberglass Translucent Conical Post with enhanced gripping functionality for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,8 to 1,4 mm. Upper: $\varnothing$ from 1,25 to 1,85 mm. Tapered and knurled shape.
TPOV 15 TPOV 17 TPOV 19	Fiber Post Glass Translucent Oval 	a 0,9 - b 1,5 a 0,9 - b 1,7 a 0,9 - b 1,9	Translucent Fiberglass Oval Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs.	L = 15 mm. a : short axis = 0,9 mm. b : long axis from 1,5 to 1,9mm. Tapered shape.
CPCO2-4 12 CPCO2-4 14 CPCO2-4 16 CPCO2-4 18 CPCO2-4 40	Fiber Post Conical Carbon 2-4% 	Conicity 2-4% $\varnothing$ 0,5 - 1,2 $\varnothing$ 0,7 - 1,4 $\varnothing$ 0,9 - 1,6 $\varnothing$ 1,1 - 1,8	Carbon Fiber Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,2 to 1,8 mm. Tapered shape 2-4%.
CPCO4-6 15 CPCO4-6 17 CPCO4-6 19 CPCO4-6 21 CPCO4-6 40	Fiber Post Conical Carbon 4-6% 	Conicity 4-6% $\varnothing$ 0,5 - 1,5 $\varnothing$ 0,7 - 1,7 $\varnothing$ 0,9 - 1,9 $\varnothing$ 1,1 - 2,1	Carbon Fiber Conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. Lower: $\varnothing$ from 0,5 to 1,1 mm. Upper: $\varnothing$ from 1,5 to 2,1 mm. Tapered shape 4-6%.
WPCORET 10 WPCORET 12 WPCORET 14 WPCORET 16 WPCORET 18 WPCORET 40	Fiber Post Conical Glass RETENTIVE 	$\varnothing$ 1,0 $\varnothing$ 1,2 $\varnothing$ 1,4 $\varnothing$ 1,6 $\varnothing$ 1,8 $\varnothing$ 1,0 $\varnothing$ 1,2 $\varnothing$ 1,4 $\varnothing$ 1,6	Fiberglass Conical Post with enhanced retentive functionality for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. $\varnothing$ from 1,0 to 1,8 mm. Tapered shape.
CPCORET 10 CPCORET 12 CPCORET 14 CPCORET 16 CPCORET 18 CPCORET 40	Fiber Post Conical Carbon RETENTIVE 	$\varnothing$ 1,0 $\varnothing$ 1,2 $\varnothing$ 1,4 $\varnothing$ 1,6 $\varnothing$ 1,8 $\varnothing$ 1,0 $\varnothing$ 1,2 $\varnothing$ 1,4 $\varnothing$ 1,6	Fiberglass Conical Post with enhanced retentive functionality for endodontic reconstruction. Radiopaque. Micro-rough surface. Supplied in boxes of 10 pcs. Kit 40 pcs. (10 for each size)	L = 21 mm. $\varnothing$ from 1,0 to 1,8 mm. Tapered shape.

* Each diameter post has the correspondent reamer

ENDODONTIC DRILL - REAMERS

CODE Price	PRODUCT	SIZE	DESCRIPTION	TECNICAL INFORMATION
CYDR 08 CYDR 10 CYDR 12 CYDR 14 CYDR 17	Cylindrical Reamer 	ø 0,8 ø 1,0 ø 1,2 ø 1,4 ø 1,7	Stainless steel drill for cylindrical post. The drill is guaranteed for more than 30 sterilization cycles without losing its cutting capabilities or becoming dark. It is designed for a lateral cut (root canal) with a profile few microns more than the referred post to allow cementation.	The cylindrical drill matches the following criteria to match the post diameter size: No ring 0,8 mm. 1 ring ø 1,0 mm. 2 rings ø 1,2 mm. 3 rings ø 1,4 mm. 4 rings ø 1,7 mm.
CODR2 08 CODR2 10 CODR2 12 CODR2 14	Conical Reamer 	Conicity 2% ø 0,8 - 1,25 ○ ø 1,0 - 1,45 ● ø 1,2 - 1,65 ● ø 1,4 - 1,85 ●	Stainless steel drill for conical post. The drill is guaranteed for more than 30 sterilization cycles without losing its cutting capabilities or becoming dark. It is designed for a lateral cut (root canal) with a profile few microns more than the referred post to allow cementation.	The conical drill matches the following: White: ø 0,8 - 1,25 mm. Yellow: ø 1,0 - 1,45 mm. Blue: ø 1,2 - 1,65 mm. Black: ø 1,4 - 1,85 mm.
CODR2-4 12 CODR2-4 14 CODR2-4 16 CODR2-4 18	Conical Reamer 	Conicity 2-4% ø 0,5 - 1,2 ○ ø 0,7 - 1,4 ● ø 0,9 - 1,6 ● ø 1,1 - 1,8 ●	Stainless steel drill for conical post. The drill is guaranteed for more than 30 sterilization cycles without losing its cutting capabilities or becoming dark. It is designed for a lateral cut (root canal) with a profile few microns more than the referred post to allow cementation.	The conical drill matches the following: White: ø 0,5 - 1,2 mm. Yellow: ø 0,7 - 1,4 mm. Blue: ø 0,9 - 1,6 mm. Black: ø 1,1 - 1,8 mm.
CODR4-6 15 CODR4-6 17 CODR4-6 19 CODR4-6 21	Conical Reamer 	Conicity 4-6% ø 0,5 - 1,5 ○○ ø 0,7 - 1,7 ●● ø 0,9 - 1,9 ●● ø 1,1 - 2,1 ●●	Stainless steel drill for conical post. The drill is guaranteed for more than 30 sterilization cycles without losing its cutting capabilities or becoming dark. It is designed for a lateral cut (root canal) with a profile few microns more than the referred post to allow cementation..	The conical drill matches the following: White: ø 0,5 - 1,5 mm. Yellow: ø 0,7 - 1,7 mm. Blue: ø 0,9 - 1,9 mm. Black: ø 1,1 - 2,1 mm.



kit with drills



Kit

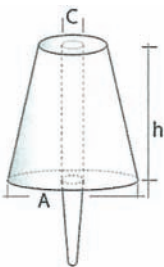


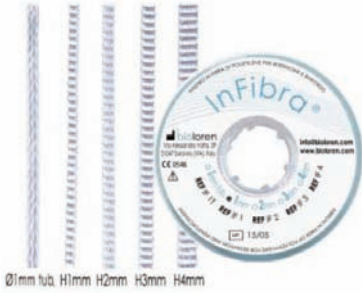
Refil



Kit drills

SPECIAL POSTS

CODE Price	PRODUCT	DESCRIPTION	TECNICAL INFORMATION
MONOCORE S MONOCORE L MONOCORE START KIT MONOCORE KIT		Monocore is a fiberglass post made of a upper part (abutment side) and a lower conical part (post side). The upper part is easily reshaped with a diamond bur for a quick reconstruction. The lower part can be inserted in a root canal or in an implant and is self-threading. Micro-rough surface. Refill: 3 small or 3 large START KIT: 3 small + 3 large KIT: 3 small + 3 large + 2 reamers	Monocore Small Conical part: $\varnothing$ from 0,7 to 1,4 mm. L post = 7 mm. Monocore Large Conical part: $\varnothing$ from 1,1 to 1,7 mm. L post = 9 mm.
ADAPT		Adapt Post is a fiberglass post made of a upper part (abutment side) and a lower conical part (post side). The upper part is easily reshaped with a diamond bur for a quick reconstruction. The lower part can be inserted in a root canal or in an implant and is self-threading. Micro-rough surface. Supplied in kit 3 pcs without drill.	Adapt Post Conical part: $\varnothing$ from 1,7 to 2,1 mm. L post = 7,0 mm. Abutment $\varnothing$ 6,0 mm.
FIT		Fit Post is a fiberglass post made of a upper part (abutment side) and a lower conical part (post side). The upper part is easily reshaped with a diamond bur for a quick reconstruction. The lower part can be inserted in a root canal or in an implant. Micro-rough surface. Supplied in kit 3 pcs without drill.	Fit Post Conical part: $\varnothing$ from 1,0 to 1,4 mm. L post = 7,0 mm. Abutment $\varnothing$ 4,5 mm.
POST BALL S POST BALL L POST BALL Kit		Post ball is a white fiberglass conical Post for endodontic reconstruction. Radiopaque. Micro-rough surface. It is provided together with a small cap to connect the ball-side. Refill 4 small or 4 large Kit: 2 small + 2 large.	Post Balll L = 21 mm. Small Ball: $\varnothing$ 1,8 mm. Large Ball: $\varnothing$ 2,5 mm. Conical shape.
Bicore		Bicore is made of two separate parts: a post part and a stump part. It offers very good aesthetic results in case of direct restorations or prosthetic crown. It has a better capacity to absorb forces. 3 stumps 3 posts	Bicore A = 6 mm. H = 8 mm. C = $\varnothing$ 1,4 mm. Post: $\varnothing$ 1,4 mm.

CODE Price	PRODUCT	SIZE	DESCRIPTION
25 cm.: IF 1T IF 1 IF 2 IF 3 IF 4 50 cm.: IF 1T IF 1 IF 2 IF 3 IF 4	InFibra Ribbon 	Ribbon for Periodontal Splinting, Spacers and Retention, Fixing of Bridges and Prostheses, Reinforcement of Loose Teeth, Orthodontic Retainers. Sold in rolls of 25 or 50 cm.	InFibra ribbon is made in polyethylene. Height: 1 mm tubular 1 mm ribbon 2 mm ribbon 3 mm ribbon 4 mm ribbon
IF 25	InFibra Large Ribbon 	Large Ribbon especially used in dental laboratories. Internal and external frameworks for bridges, Reinforcement of full or partial dentures, Fixing of Bridges and Prostheses. Sold in rolls of 50 cm.	InFibra large ribbon is made in polyethylene. Height 25 mm.
IF KIT 5	InFibra Ribbon Kit 5 	Ribbon for Periodontal Splinting, Spacers and Retention, Fixing of Bridges and Prostheses, Reinforcement of Loose Teeth, Orthodontic Retainers. Sold in roll of 50 cm. Provided with or without special scissors for net cut. Packaging in plastic boxes, supplied with 5 rolls.	Height: Ø 1mm.Tubular H. 1 mm ribbon H 2 mm ribbon H 3 mm ribbon H 4 mm ribbon
SCS	Professional Scissor 	Scissors: Special scissors for net cut. Packaging: box 1 pcs	Features: Micro-shagreened for a net cut.



InFibra Ribbon Kit 5 rolls



**InFibra
H 25 mm. rolls 50 cm.**



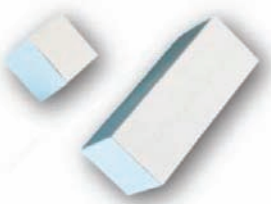
**Ribbon for retention
and reinforcement**

COMPOSITES

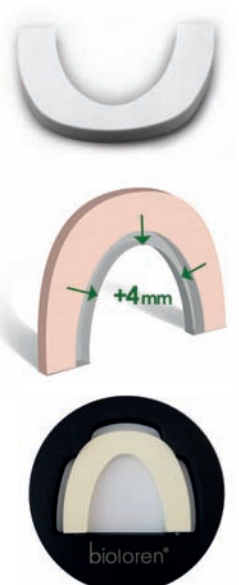
CODE Price	PRODUCT	SIZE	DESCRIPTION
Universal Resin Cement URC 7 L URC 7 M		U.R.C. Universal Resin Cement. It is resin based cement is ideal for cementation of fiber posts, crowns, bridges and inlay, metal, zirconia and Polyethylene. Dual and hydrophilic. It forms a positive associations with water, taking advantage of the moisture that is always present in the mouth. Packing: automix syringe 7 grams 7 gr syringes + 10 straight or curved mixing tips:	Self-adhesive: no primers, no silane, no bonding agents required. Self etching to dentin. Colour stable when cured. Available in: Medium viscosity. Low viscosity When activated the cured material is no more affected by water. Compressive Strength: 310 Mpa Tensile Strength: 55 Mpa Retention value: 32 kg.
URCP D URCP A		Straight mixing tips. Angled mixing tips.	Packing: 10 pcs.
KIT INTRO ARCH		Kit to make complete reconstruction of 2 maxillaries. It is used manually only by a micromotor	Packing 1 urc 7 g + 5 mixing tips 2 arches with different size
LC Bonding IBLTP		IBLTP 5 ml bottle, 5 pcs micro applicators, cup. I-Bonding LC is a 5 th generation (il th mettito in alto io non sono capace) generation light curing, bioinert, single component one-step adhesive, suitable for use with wet bonding technique. Etching of the surfaces is required.	User to bond dentin and enamel to polymer based restorative materials, compomers, ceramic surfaces to improve retention of a restoration.



TRILOR® FIBER DISKS AND BLOCKS

CODE Price	PRODUCT	SIZE	DESCRIPTION
FDS 10 FDS 12 FDS 14 FDS 16 FDS 18 FDS 20 FDS 250 FDS 14 PK FDS 16 PK FDS 20 PK		TRILOR Disk is an element for CAD/CAM system. A fiberglass disk with the exclusive interlacing of the fibers and resin matrix, that provides high performances Trilor disk us the ideal material to build meso-structures and full arches. Permanent Supplied in box 1 pcs TRILOR Disk Pink (like Trilor Disk, but in pink color) Trilor disk pink allows to solve some prosthetic or clinical needs	ø 98 mm. H 10 mm. H 12 mm. H 14 mm. H 16 mm. H 18 mm. H 20 mm. H 25 mm. H 14 mm. H 16 mm. H 20 mm.
FDZ 14 FDZ 16 FDZ 18 FDZ 20 FDZ 25 FDZ 14 PK FDZ 16 PK FDZ 20 PK		TRILOR Disk is an element of the CAD / CAM system. It is a fiberglass disc with a special weave of fibers and epoxy resin that offers high performance. Tread disc is the ideal material for constructing substructures and complete arches. Permanent. Supplied in single plastic box. Zirkonzahn compatible *. TRILOR Disk Pink (like Trilor Disk, but in pink color) Trilor disk pink allows to solve some prosthetic or clinical needs	ø 95 mm. H 14 mm. H 16 mm. H 18 mm. H 20 mm. H 25 mm. H 16 mm. H 20 mm.
FBS 201514 FBS 201915 FBS 401915 FBS 551915 FBS 652522 FBS 654022 FBS 854022		TRILOR Block is an element of the CAD / CAM system. It is a glass fiber block with a special weave and epoxy resin that offers high performance. Permanent. Supplied in single box.	mm. 20x15x14 20x19x15 40x19x15 55x19x15 65x25x22 65x40x22 85x40x22

TRILOR® ARCH

CODE Price	PRODUCT	SIZE	DESCRIPTION
TA 3.5 TA 5.5 TA 7.5 TA KIT (1 for each size) TA 3.5 NW PK TA 5.5 NW PK TA 7.5 NW		TRILOR ARCH is made of glass and resin fibers. Indications for use: permanent reinforcement bar, for immediate loading, implant position transfer key, implant surgical guide Ideal for All On four and All on Six, or for small jobs such as emiarccates or bridges. Permanent. Supplied in box 3 pcs MANUAL USE, ONLY MICROMOTOR NO CAD/CAM MILLING MACHINE TRILOR ARCH Starter Kit was created to introduce clinicians and technicians to the use of the Trilor Arch Metal Free techno-polymer. The adapter disc allows the Trilor Arch 7.5 mm arch to be positioned in order to work it with CAD/CAM milling machines. The adapter disc will be inserted into the CAD/CAM machine.	H 3,5 mm. H 5,5 mm. H 7,5 mm. H 3,5 mm. H 5,5 mm. H 7,5 mm. 1 Adapter + 2 Arch

* Zirkonzhan is a trademark of Zirkonzhan spa

ON REQUEST - PPE Personal Protection Equipment



SURGICAL MASK 50 pcs
€ 3,00

3 Layers
Water resistant
Anatomic design
Filtering efficacy >99%
Box 50 pcs



FFP2 MASK 25 pcs
€ 6,00

Filter mask FFP2 NTI
Single packed
Box 25 pcs
Anatomic design



BIS - Bioloren Immune Shield
€ 7,00

Protective vision anti-droplet
Scratch proof in PTE
Single packed

Trilor[®]

Arch

For dental use only

3 PCS


☐ **3,5mm**
 REF TA 3,5

☐ **5,5mm**
 REF TA 5,5

☐ **7,5mm**
 REF TA 7,5

☐ **KIT**
 REF TA KIT

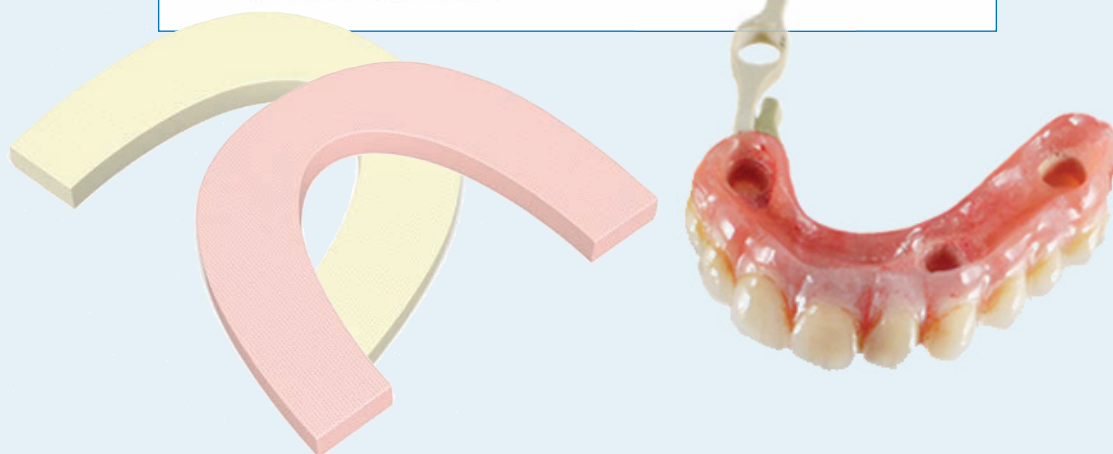

Via Alessandro Volta, 59
 21047 Saronno (Va) ITALY
 Ph./Fax. +39 0296703261
www.bioloren.com



0546

LOT

Rev.3



bioloren[®]
 metal free dental solutions



Bioloren Srl - Via Alessandro Volta, 59 - 21047 Saronno (VA) - Italy
 Ph./Fax: +39 02 96703261 - info@bioloren.com - www.bioloren.com

MADE IN ITALY

